



DIGITAL CAMERA ADAPTOR

TD-2



TD-2 does not include a camera.

Light path changeable in accordance with right or left better eye through one-touch operation

The light path to capture images into the digital camera can be changed in accordance with the examiner's better eye. By adjusting the dial on the back, you can change the light path through one-touch operation. The fundus can be observed without worrying about displacement between an image in the field of view and a captured image.



Yellow filter IN/OUT mechanism

High-contrast images can be obtained by using a built-in yellow filter in cornea fluorescence observation. Easy filter IN or OUT is available through one-touch operation. In its setting with our slitlamp, high-contrast images can be obtained by the best matching with a cobalt blue filter.

Camera models compatible with TD-2

TD-2 types do not come with a digital camera. For the purchase of a digital camera, refer to the following compatible cameras:

Compatible cameras manufactured by Canon
 EOS kiss Digital N
 EOS DIGITAL REBEL XT/350D DIGITAL
 EOS DIGITAL REBEL XT/300D DIGITAL
 EOS 10D, EOS 20D, EOS D30, EOS D60



Slitlamp models compatible with TD-2 types



Type A

Compatible models
 SM-70N, SM-70,
 SM-90N, SM-90



Type B


Compatible models
 SM-2N, SM-2,
 SM-4N, SM-4

Capturing of still image with one push

By connecting a trigger cable between the slitlamp and camera, a still image can be captured every time the joystick button of the slitlamp is pushed.



The digital camera adaptor TD-2 allows you to import the images taken into the PC and display them on the monitor by installing a commercially available single-lens reflex digital camera and combining it with a slitlamp.

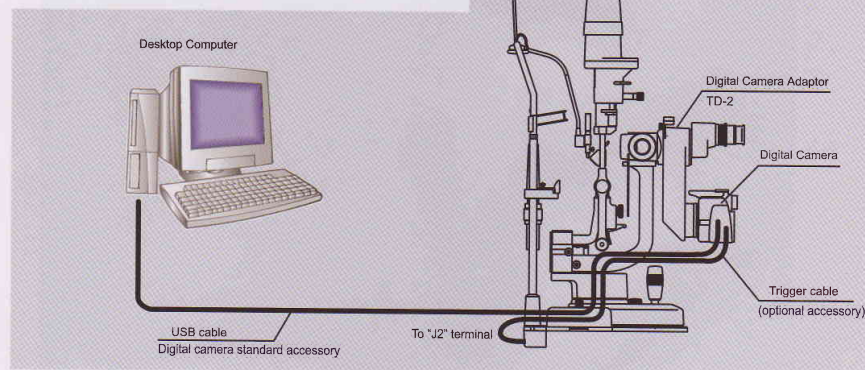

Optimed Pty Ltd

Unit 1, 15-17 Chaplin Drive, Lane Cove, 2066
 Ph: 1300 657 720 Fax: 02-9420 1144
 sales@optimed.com.au

iming at new levels in quality

TAKAGI SEIKO CO.,LTD.

330-2 IWAFUNE, NAKANO-SHI, NAGANO-KEN, 383-8585, JAPAN
 TEL.+81-269-22-4512 FAX.+81-269-26-6321
 URL: http://www.takagi-j.com E-mail: info@takagi-j.com



• Design and specifications are subject to change as improvements are made to the product.



Our quality system is certified for ISO 13485.
 R08002 Rev.0 Printed in Japan 4 2002-10-1